


Waysafe SA640

Sample Analysis Enclosure

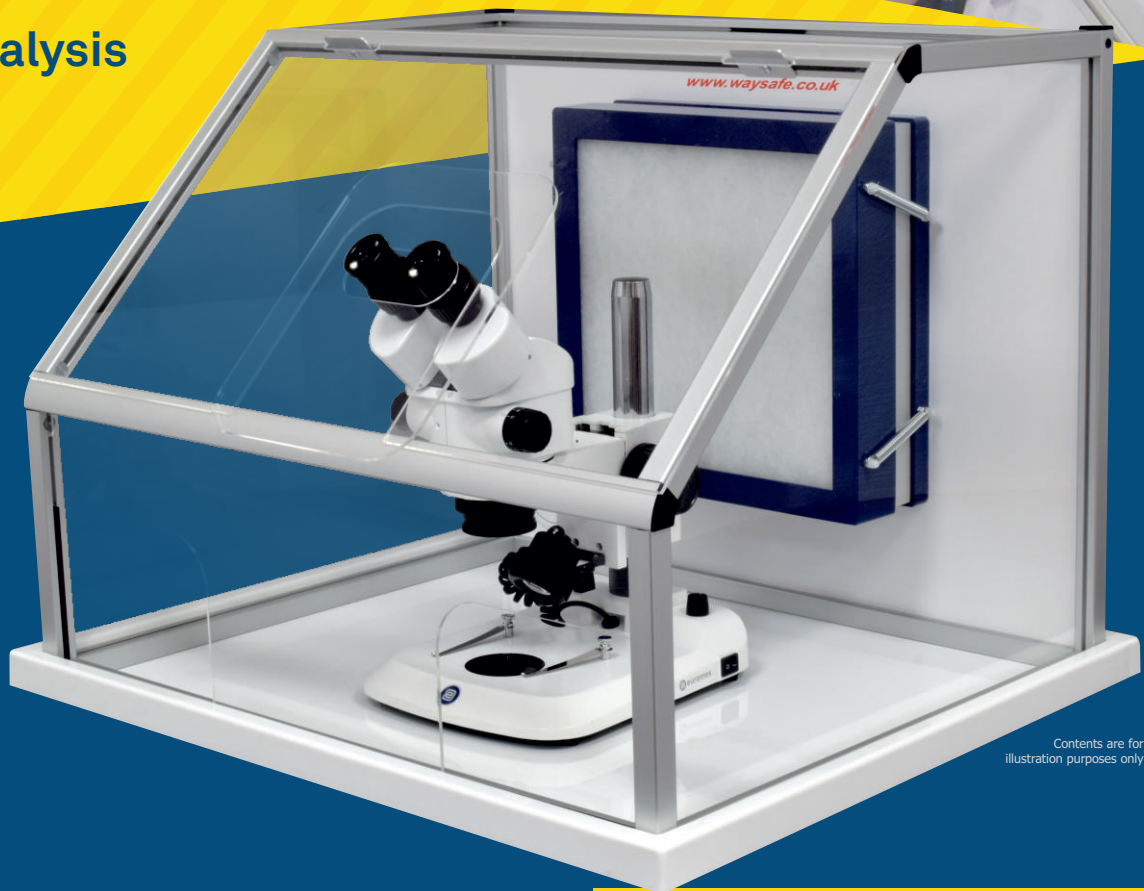

Contents are for illustration purposes only

At 640mm wide the SA640 (previously Waysafe 2) has ample space for most types of microscopes, yet is still small enough to fit on a standard benchtop.

A low noise 2 speed fan draws air through our trusted 300 series H14 HEPA filter before returning it to the lab.

The front panel can be easily removed to provide access to the interior for the placing and positioning of equipment. The lower panel has two hand holes for operator access.

Various options available including a range of filters, 120VAC option and an integral composite aluminium base in a range of colours.



- Ductless - air is filtered and clean air circulated back into the workplace
- H14 HEPA filtration (99.999% efficient)
- Easy benchtop installation
- Low maintenance
- 2 speed fan
- Spill tray provided as standard
- Lift off front for ease of access
- Fast global delivery
- Protects operators
- Perfect for microscopy and all types of sample analysis including Asbestos
- Can also be used for other processes media weighing etc.
- Made to order at our factory in the UK

Call us on +44(0)161 339 6187 or email info@waysafe.co.uk for more information

www.waysafe.co.uk

Technical Specification

Waysafe SA640

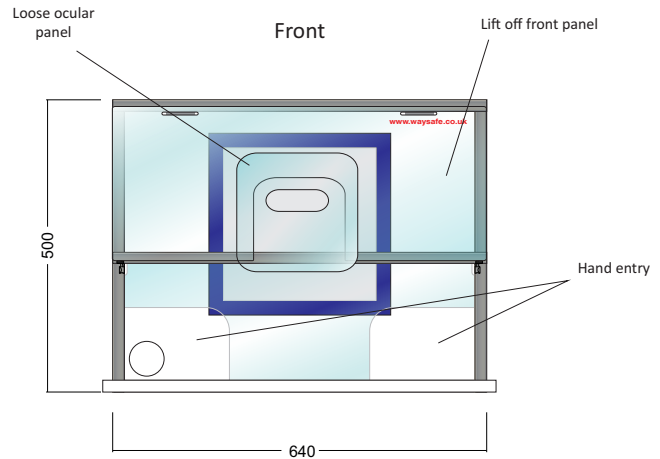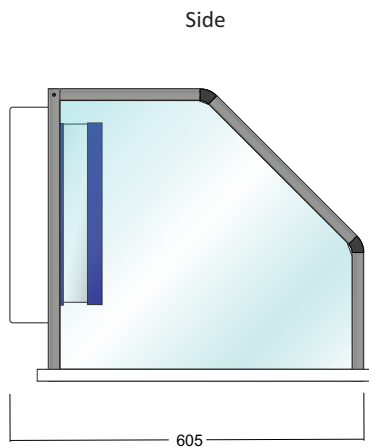

Technical Specification

Electrical	230VAC 50w, meets LVD and is CE Marked (120VAC available)
Air Speed	High > 0.5m/s, Low 0.3m/s (Face velocity)
External Dim.	Height 500mm, depth 605mm, width 640mm (Tray dimensions 640mm x 680mm)
Filter	Waysafe '300 series' H14 HEPA filter, 99.999% efficiency @ 0.3 micron
Pre filter	Fitted as standard (simple low cost G4 grade filter pad protects the HEPA filter from large particles to extend service life)
Base tray	Push fit moulded spill tray (see options)
Weight	13.9Kg with standard HEPA filter

Material Specification

Base	High Impact Polystyrene
Glazing	Cast Acrylic Glazing
Frame	Durable Anodised Aluminium
Filter Frame	Powder Coated Aluminium Frame
Ancillaries	Rubber Seals, Durable Hinges

Optional Extras

Built-in Base	Various colour options including blue, red, white and stainless steel
Left or Right Side Door	Full size opening side door (Only available with Built-in base)
Face Panel	Bespoke panels available to suit varying equipment types
Light	Waysafe TL640 lightbox
Filters	Various filter options include HEPA for particulates and carbons for odours
Variable Speed Fan	Rotary control allows users to steplessly adjust fanspeed for optimum airflow
Hours Run Meter	Electro-mechanical, non-resettable run meter for filter and service monitoring



Call us on +44 (0)161 339 6187 or email info@waysafe.co.uk for more information

www.waysafe.co.uk

